

## Technical Service Bulletin # 1

**Product:** Incubator 8000 SC


**Date:** November 30, 1999

**Re:** Incubator 8000 SC, increase of operational security of the fan

**Reference Doc:** IDM no. 1 / 06.91

**Reason:** Fan failure alarm has occurred sporadically in some Incubator 8000 SC devices.

**Solution:** In order to improve the fan monitoring in the Incubator 8000 SC, a field change must be implemented as a precautionary measure. A screw (M3) is used both as coil core and mount. Three washers must be used to achieve the proper distance between the screw head and fan wheel.

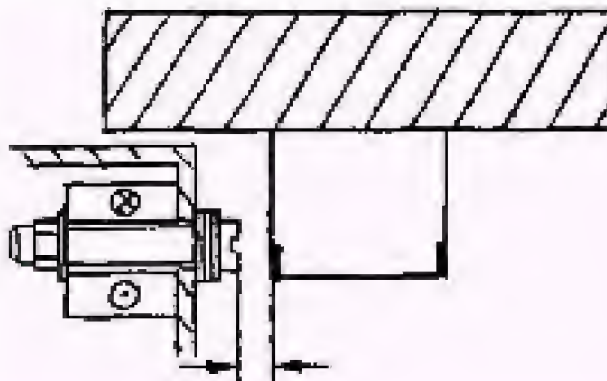

**Devices affected:** Incubator 8000 SC shipped prior to January 1991 with the following serial numbers:  
ARCM-0001 to ARCM-0015, and  
ARCN-0001 to ARCN-0015

**When:** At next Service Call or Preventive Maintenance

**Additional Info:** N/A

**Cost:** Covered by Customer, post warranty.

**Ordering Info:** possibly 3 washers A3.2 DIN 125 with the P/N 13 30 233.

**Distribution:** Dräger Service Personnel and Authorized Service Organizations for NCS products.

If you have any questions, please contact Technical Support by phone at 1-800-543-5047 or by fax at 1-215-723-5935

Dräger Medical, Inc.  
Technical Product Manager

Andreas Lenke